

USER MANUAL

ABOUT BOOSTER PRO X2

PRODUCT OPERATION

01. Press and hold the ③ Power Button for 2-3 seconds to turn on the device. When powered on, the ② LCD Screen lights up with the ON sound, and the device automatically starts in Booster X2 Mode at Level 1.

02. Briefly press the ④ Mode Button once to select the desired mode. Each press changes the mode, which is shown on the ② LCD Screen along with a voice guide and the mode-specific LED color on the ⑤ Central LED Indicator.

*Modes change in the following order: Booster X2 Mode → MC X2 Mode → Derma Shot X2 Mode → Dual Mode → Mask Mode → Air Shot X2 Mode.

03. After selecting the mode, briefly press the ③ Level Button once to adjust the level. Each press increases the level by one, and the changed level is displayed on the ② LCD Screen each time the level increases.

*You can set the level from 1 to 6, and it is recommended to start with level 1 when using the device for the first time.

04. After selecting the mode and level, use the device according to the usage method for each mode. Be sure to familiarize yourself with the usage method for each mode before use.


05. During use, the usage time is announced by voice guide every 1 minute. When the usage time reaches 5 minutes (10 minutes for Mask Mode), a care completion voice guide will play, and the mode standby screen will appear on the ② LCD Screen for 10 seconds. If there is no additional button/app operation within the 10-second countdown, the device will automatically turn off. If you would like to continue with the next mode, briefly press the ④ Mode Button once within 10 seconds. If you would like to just continue the care, briefly press the ③ Power Button once within 10 seconds.

06. Press and hold the ③ Power Button and the ④ Mode Button together for more than 3 seconds to change the voice guide language. English, Korean, Japanese, and Chinese are available.

07. For hygienic use, be sure to clean the ① Head with a dry cloth or cotton pad after use and store it properly.

08. To charge the device, connect the provided Type-C charger to the ⑦ Charging Port on the lower rear side of the device.

*When charging, a charging connection sound will play, and the charging status will be displayed on the ② LCD Screen. The device will not operate while charging.

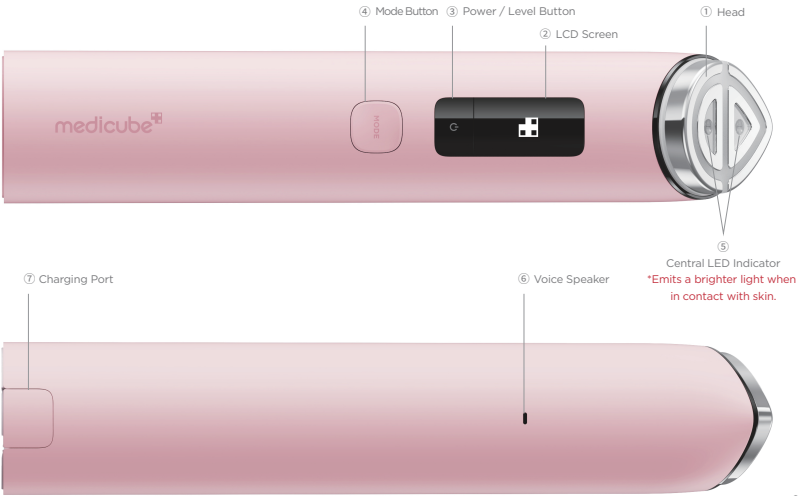

01

PRECAUTIONS FOR BLUETOOTH CONNECTION

01. Remote Control
The remote control functions operate only when the device is connected. LED color change, volume control, vibration ON/OFF, Dual Mode change, and AI Custom Care functions are available only within the remote control after Bluetooth connection. The volume level, vibration ON/OFF, Dual Mode setting, and LED color set in the remote control are saved.

02. Automatic Report
Usage records created while the device is used without being connected to the app are temporarily stored in the device. When the device is connected to the app, the data is automatically uploaded to the app report. The usage level and LED color for each mode in the report are saved as the longest-used level and color.

03. AI Care
Personalized care is provided based on the user's usage pattern and skin needs. AI Custom Care is available only when connected to the app.

04. Default LED Color Changes for Each Mode
The LED color for each mode can be changed through the remote control only when the device is connected to the app. The last used LED color is saved as the default color.

05. Device Connection Failure
Ensure that the device is not connected to another smartphone. If it is connected to another smartphone, disconnect it from the connected smartphone's app by going to [Connect Booster Pro X2] → [Settings] → [Device Disconnect], then try again.

06. Changing the Connected Smartphone
Run the app on the currently connected smartphone and disconnect the device through [Connect Booster Pro X2] → [Settings] → [Device Disconnect], then try connecting with a new smartphone.

Available only with Bluetooth connection

• Mode-Specific Color Change

• Change Dual Mode

• Sonic Vibration ON/OFF

• Volume Control

• Check Care Report

• AI Care

PRECAUTIONS

Precautions during Use

01.Be careful not to drop or shock the device, as it may cause malfunction or damage.

02.Ensure that no water or cosmetics enter the inside of the main unit.

03.If the product is not operating correctly, do not disassemble it; instead, contact customer support.

04.Remove all metal accessories in contact with the user's body.

05.Do not use the device in conjunction with electrical appliances such as electric heating pads.

06.Do not operate the device in hot places such as near fire or direct sunlight, and do not place the main unit or cable near heating appliances.

07.Do not keep it fixed on one area for an extended period during device use; make sure to move it continuously.

08.Contact with other products may cause malfunction. Please use this device with caution.

09.Do not use the device for purposes other than its original intention; it may lead to accidents.

10.Do not use this device if you regularly suffer from paralysis or have difficulties expressing yourself.

11.Do not tap the head that contacts the skin directly on hard surfaces; it may lead to accidents.

12.Do not strike or drop the main unit on hard objects.

13.Avoid using the device in high temperature, high humidity, or dusty environments, such as bathrooms, kitchens, and hot cars.

14.When using, maintain a distance from electrical appliances such as other beauty devices, TVs, and radios.

15.Do not use the device on eyelids, Adam's apple area, wounds or areas where procedures or surgeries have been performed.

16.Do not use the device at voltages other than the rated voltage; use of improper voltages can cause heat generation, leading to fires, electric shock, or malfunction.

17.Do not handle the device or detach the charging cable with wet hands; it can lead to electric shock.

18.Do not attempt to repair, disassemble, or modify the main unit yourself; it can cause serious accidents.

19.The device will not operate while charging.

20.For accurate skin recognition, please grasp the main part of the device as if surrounding it.

21.Children should be supervised to ensure that they do not play with the appliance.

22.The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction.

Storage Precautions

01. Store in a place out of the reach of children and away from moisture.

02. Do not store in places where magnetic fields or electromagnetic waves are generated; they may cause malfunction.

03. After use, clean the head with a wet wipe, clean cloth, or soft towel and store it dry.

04 Ensure that the device's power is turned off before storage.

05. Do not use flammable substances such as paint thinner, gasoline or benzene for cleaning, as it can cause chemical reactions, discoloration, and deformation.

User Precautions

Do not use if you are in poor health, under the influence of alcohol, taking sleeping pills, experiencing fatigue, or are under medical treatment. If you experience pain or unusual symptoms during use, discontinue use immediately. **Repeated use of the device on the same area or improper usage may cause damage to the skin. If any abnormality appears on the skin, cool the affected area and consult a medical professional immediately.**

* Even in regular use, individuals with sensitive skin may experience slight redness. Please read the instructions and precautions carefully and perform proper aftercare.

Do not use if you have any of the following symptoms or if you fall under any of the following conditions.

01. Individuals with symptoms of redness and stinging due to excessive sun tanning 02. Individuals with allergies, especially those with sensitive skin 03. Individuals using electronic devices worn directly inside the body, such as cardiac pacemakers 04. Individuals using electronic medical devices for life support, such as artificial heart-lung machines 05. Individuals using wearable medical devices such as electrocardiogram devices 06. Individuals taking steroid-based hormone medications for an extended period 07. Individuals with liver function disorders who have symptoms of telangiectasia 08. Individuals undergoing outpatient treatment 09. Areas with inflamed acne 10. Areas with dermatitis caused by cosmetics, or individuals with atopic dermatitis 11. Areas with blemishes accompanied by itching, or areas

02

with blemishes caused by physical irritation 12. Areas with metal, plastic, silicone, or other foreign objects inside the body 13. Patients with blood disorders, heart diseases, acute illnesses, bleeding areas, malignant tumors, fractures, infectious diseases, eyelids, trauma areas, skull areas 14. Children, pregnant women, breastfeeding women, and elderly individuals 15. Individuals who regularly suffer from paralysis or have difficulties expressing themselves. * This product is intended for cosmetic use only. It is not intended to diagnose, treat, cure or prevent any medical condition or disease. Please read all instructions and warnings before use. Use only as directed. Always consult with a healthcare professional if you have any medical concerns or conditions before using the device.

SPECIFICATION

Model Name ME-BPXT-K **Product Name** MEDICUBE AGE-R BOOSTER PRO X2
Seller APR US INC **Contact** medicube_uscs@apr-in.com **Seller Address** 41 GREENFIELD IRVINE CA 92614
Manufacturer APR Factory Co., Ltd.
Manufacturer Address 2F, 107 Gasan Digital 2-ro, Geumcheon-gu, Seoul, Republic of Korea/ 50 Wonseokjeongyongso-gil, Poseung-eup, Pyeongtaek-si, Gyeonggi-do, Republic of Korea
KC Registration No. R-R-pg1-ME-BPXT **FCC ID** 2AC7Z-ESPS3WROOM1
Instructions and Precautions Refer to the User Manual
Serial Number / Date of Manufacture Refer to the Specification Label
Other Required Information Refer to the User Manual

Product Specifications
Dimensions 39 x 195 mm (1.53 x 7.68 inches) **Input Voltage** DC 5V / 2A (USB-C) **Battery** 3.7V / 2,600mAh Li-ion **Power Consumption** max. 3.5 Watt (3.7V x 0.95A) **Weight** 230g (0.51lbs)
Charging Time Approx. 4 hours
MADE IN KOREA
Intended Use Facial Beauty Massaging Device

Battery Specifications
Quality Indication in accordance with the Electrical Appliances and Consumer Products Safety Control Act
Safety Confirmation Report No.: ZUI0601-25001

K **Model Name** ND18650-2600mAh-3C **Product Name** Battery (LI-ion Battery)
Rating 3.7V / 2,550mAh **Manufacturer** DONGGUAN LEAN POWER NEW ENERGY TECHNOLOGY CO.,LTD **MADE IN CHINA**

This device complies with Part 15 of the FCC rules.
Operation is subject to the following two conditions:
(1) This device may not cause harmful interference, and (2) this device must accept any

interference received, including interference that may cause undesired operation.
This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:
- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



FCC Note Caution : Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

EXCHANGE / RETURN POLICY

Warranty
For all MEDICUBE devices, we offer a 1-year limited warranty period from the date of your original purchase. If your device does not work properly due to a manufacturer's defect, we will replace the defective device free of charge within 1 year. This warranty does not cover the devices damaged by: 01. Accident, misuse, or alteration 02. Use with unauthorized accessories or use other than instructed 03. Connecting to incorrect current and voltage. This warranty only extends if the devices are used in conjunction with authorized device accessories (such as the enclosed charging cables).
Returns
For the MEDICUBE devices, all returns must be initiated within 30 days of receiving your order. Please refer to the information as shown below. **Simple Remorse:** A refund due to a change of mind is only applicable within 30 days of receiving and will only be accepted if the product is returned in original, sealed, and unopened condition. *In the case of a return for simple remorse (change of mind), the customer needs to bear the shipping fee. For international orders: You will be responsible for paying for customs fee and returning shipping fee. Additionally, we cannot provide free return shipping labels for international customers. **Damaged or Defective Product:** 1) Please reach out to us via our Customer Service Center email with the following details: First name and last name, email address, phone number, order number. Kindly attach a photo or video describing your issue. 2) We will proceed to issue a confirmation email for you (please note that our Customer Service usually takes 2-3 business days to respond). 3) After receiving a confirmation email, please send the product back to our warehouse via the return address. You will be informed of the return address after the email confirmation. 4) Once the product has arrived at our warehouse, our staff will inspect the condition of the product. 5) In case of product defect, we will proceed to exchange and send a new product to you. *Note: This process may take up to 1 month, including the time taken to ship the product back to our warehouse, and the time taken to inspect the product.

Subject of Recommendation of Lithium Battery Discarding and Charging Method
Warning: Since it may accompany danger such as heating, fire, and explosion, please observe the following matters. **1) Handling Matters** A visibly swollen battery may be dangerous, immediately contact the manufacturer or the seller. Use only the designated charger. Keep away from fire (Do not put inside a microwave oven). Do not leave it inside a car in summertime. Do not store or use at high temperature and humid place such as sauna. Do not use it for a long time while placing it on blanket, electric pad, and carpet. Do not store at a sealed space while leaving the power on. Be cautious not to contact the electrical terminal with metal products such as necklace, coins, key, and watch. Use the secondary lithium battery warranted by the manufacturer. Do not implement behaviors such as disassembly, compressing, and penetrating. Do not apply abnormal impact such as dropping it from a high place. Do not expose to a temperature above 50 degrees Celsius. Do not contact with humidity. **2) Disposal Guideline** Do not dispose of used batteries carelessly. When disposing the battery, do not incinerate or throw batteries into fire. Heat may cause explosion or fire. The method of battery disposal may vary by country and region. Dispose with appropriate method. **3)Charging Recommendations** Do not charge the battery with the method other than the one shown in the user manual. Fire or explosion may occur.
Complete Charge Time : About 4 hours
Charge time may vary depending on the battery remnant amount of portable battery or surrounding temperature or charging environment. It may change without notice for performance improvement of specifications and exterior appearance.

Please be sure to read the instruction manual and precautions before use.
Images in the manual are for illustrative purposes only.

Copyright © 2026, APR Corp. All Rights Reserved.
medicube.us